
$Flex_{\scriptscriptstyle{\mathrm{TM}}}$ interproximal finishing and contouring strips

The Patented, Safe, Flexible Design of ET Flex Allows You to Adapt the Strip to Meet All of Your Interproximal Finishing and Contouring Needs.

- ✓ Orthodontics: IPR/Tooth Slenderizing
- ✓ Composite Restorations: Contouring/Finishing all Interproximal Surfaces
- ✓ Crown & Bridge Restorations: Cement Removal and Clean-Up Following Cementation



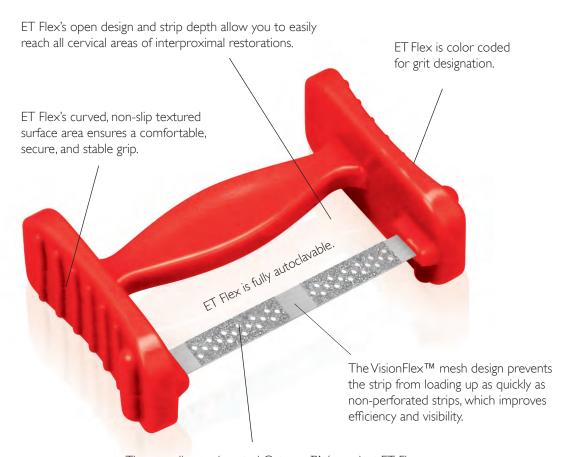

Grip the ET Flex at or above the crossbar and the strip becomes taut for contact areas.




Grip the ET Flex below the crossbar and the strip curves allowing you to perfectly adapt to any cervical or contact area contour.

The non-diamond coated Gateway™ (on select ET Flex strips) allows you to safely pass through tight contacts for interproximal finishing of composite restorations and cement clean up in the cervical area.

By Your Side in Dentistry



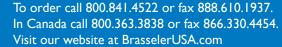





ET FLEX™ SINGLE-SIDED RESTORATIVE STRIPS WITH NON-DIAMOND COATED GATEWAY™

Ideal for safely passing through contacts for interproximal finishing of composite restorations and cement clean up in the cervical area. Narrow (2.5mm), 10 Per Pack

ET25XFSS

Yellow: Extra-Fine Grit (15 micron) Thickness: .08mm



ET25FSS

Red: Fine Grit (30 micron) Thickness: .10mm



ET25MSS

Blue: Medium Grit (45 micron) Thickness: .13mm



ET FLEX™ SINGLE-SIDED IPR AND GENERAL RESTORATIVE STRIPS

Ideal for safely and effectively finishing and contouring composite restorations and single tooth slenderization. Wide (3.75mm), 10 Per Pack













ET37CSS

Green: Coarse Grit (90 micron) Thickness: .18mm



ET37XCSS Gray: Extra-Coarse Grit (120 micron) Thickness: .22mm



ET FLEX™ DOUBLE-SIDED IPR AND GENERAL RESTORATIVE STRIPS:

Ideal for effecient IPR and finishing and countouring composite restorations. Wide (3.75mm), 10 Per Pack









ET37XCDS

Thickness: .40mm

ET FLEXTM SERRATED STRIP:

Ideal for working through tight contacts. Wide (3.75mm), 10 Per Pack

ET37SER

Black

No Diamond Coating Thickness: .05mm



Green: Coarse Grit (90 micron)



Gray: Extra-Coarse Grit (120 micron)





To order call 800.841.4522 or fax 888.610.1937. In Canada call 800.363.3838 or fax 866.330.4454. Visit our website at BrasselerUSA.com

ET FLEX™ ASSORTMENT KIT

Ideal tray setup! Go straight from the autoclave to the operatory! Best value!

K0223

- One (1) ETASST (includes one of each ET Flex)
- One (1) A625BLUE autoclavable organizational block

By Your Side in Dentistry



©2012 Brasseler USA. All rights reserved. B-3461-02.12-15M-RRD