


---

# RESTORATION REMOVAL

Amalgam, PFM, All Ceramic,  
Zirconia & Lithium Disilicate

---

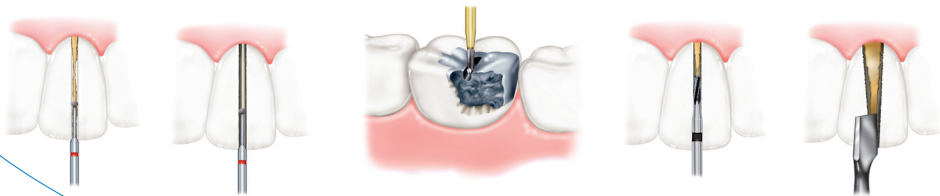


By Your Side in Dentistry



**BRASSELER USA**  
DENTAL INSTRUMENTATION

## ALLOY REMOVAL (AMALGAM)

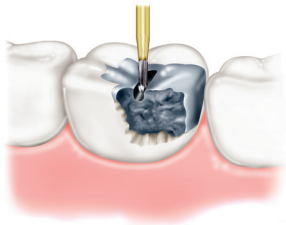

# SABRECUT



Order #	330SC	557SC	1558SC
Head Size (1/10 mm)	008	010	012


Super premium multi-purpose carbides  
in a single piece construction.

## ZIRCONIA, LITHIUM DISILICATE, ALL CERAMIC

**Dialite**<sup>®</sup>  
DIAMONDS




Order #	8801LDF	8369DF	8856DF
Head Size (1/10 mm)	014	025	016
Head Length (mm)	—	5.5	8.0

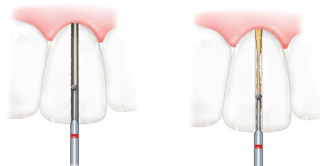

Dialite Diamonds special hybrid construction are  
ideal for removal of Zirconia & Lithium Disilicate.  
Recommended Speeds: 80-120k rpms

## METAL CUTTERS (NON-PRECIOUS)

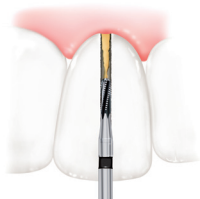

Order #	1958	H34	H34	H34L
Head Size (1/10 mm)	012	010	012	012

Ideal for removal of PFM crowns.

Tip: PFM crown removal use Dialite Diamonds to cut through ceramic, then use Metal Cutters to cut through metal substructure.

## EB134 CROWN SPREADER

Place the EB134 into the separation groove and rotate to spread and release the metal substructure from the tooth.



# Brasseler's Complete Restoration Removal System



To order call 800.841.4522 or fax 888.610.1937.  
In Canada call 800.363.3838 or fax 866.330.4454.  
Visit our website at [BrasselerUSA.com](http://BrasselerUSA.com)

**K0196**  
Zirconia/PFM Crown  
Removal System