


ANOTHER INNOVATION FROM BRASSELER USA...

VISIONFLEX™

DIAMOND/SERRATED SAW COMBO DIAMOND STRIPS

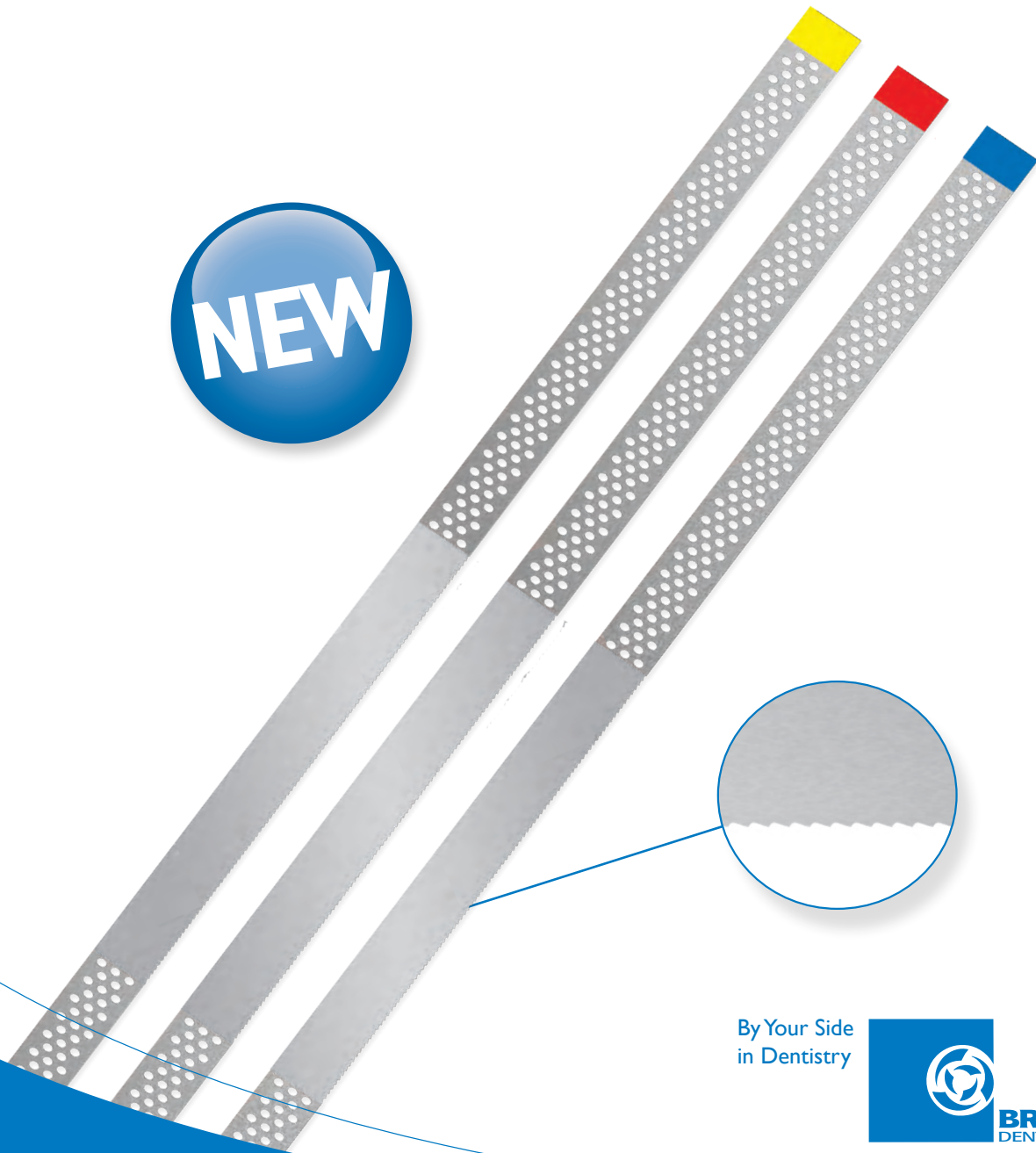

By Your Side
in Dentistry



BRASSELER USA®
DENTAL INSTRUMENTATION



ANOTHER INNOVATION FROM BRASSELER USA...

VISIONFLEX™

Diamond/Serrated Saw Combo Diamond Strips

In addition to our extensive interproximal strip line, Brasseler USA is introducing a new combination diamond and saw blade strip. This versatile strip will help meet all of your interproximal finishing and stripping needs by allowing you to move quickly from one application to the next, eliminating the need for two different strips. The serrated saw blade is ideal for opening tight interproximal contacts or cement removal, while the perforated diamond portion quickly strips enamel during orthodontics or finishes composite and porcelain restorations. Available in three grits (extra-fine, fine and medium), these strips will help you save time and money.



Order No.	Thickness (mm)	Width (mm)	Length (mm)	Saw Area (mm)	Diamond Area (mm)	Grit	Single/Double Sided	Qty/Pack
WS37EFCOMBO	0.08	3.75	147	40	53.5, 53.5	● Extra-Fine	Single	10
WS37FCOMBO	0.10	3.75	147	40	53.5, 53.5	● Fine	Single	10
WS37COMBO	0.13	3.75	147	40	53.5, 53.5	● Medium	Single	10
WS37ACOMBO	-	3.75	147	40	53.5, 53.5	Asst. (EF, F, M)	Single	15 (5 ea)

Autoclavable Diamond Strip Holder

The A101 is a fully sterilizable aluminum tube that protects strips during sterilization.



A101

Holds	Diamond Strips
Blue	A101B



To order call 800.841.4522 or fax 888.610.1937.
In Canada call 800.363.3838 or fax 866.330.4454.
Visit our website at BrasselerUSA.com

By Your Side
in Dentistry



BRASSELER USA
DENTAL INSTRUMENTATION