


solutions™

Dialite® Adjusting, Finishing
& Polishing System



Order Number:
225.00.638

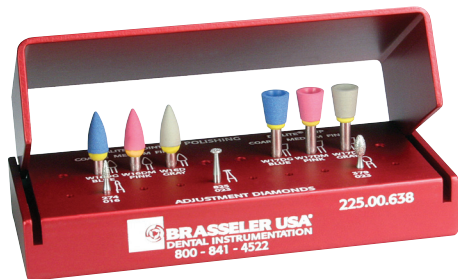





| Figure Number | WI6DG | WI6DM | WI6D | WI7DG | WI7DM | WI7D |
|--------------------------------|-------|-------|------|-------|-------|------|
| Head Size ($\frac{1}{10}$ mm) | 040 | 040 | 040 | 060 | 060 | 060 |
| Head Length (mm) | 10.0 | 10.0 | 10.0 | 7.5 | 7.5 | 7.5 |








| Figure Number | 8801LDF | 8369DF | 8856DF |
|--------------------------------|---------|--------|--------|
| Head Size ($\frac{1}{10}$ mm) | 014 | 025 | 016 |
| Head Length (mm) | — | 5.5 | 8.0 |

Only Available Direct from your Brasseler Representative
To order call 800.841.4522 or fax 888.610.1937
In Canada call 800.363.3838 or fax 866.330.4454
Visit our website at BrasselerUSA.com

solutions™ 225.00.638 Dialite® Adjusting, Finishing & Polishing System

| | Order # | Head Size | Grit | Usage |
|--|------------|-----------|----------|----------------------------------|
|  | W16DCHP-RA | 150 | C | Polisher - Occlusal Surface |
|  | W16DMHP-RA | 050 | M | Polisher - Occlusal Surface |
|  | W16DHP-RA | 030 | F | Polisher - Occlusal Surface |
|  | W17DCHP-RA | 150 | C | Polisher - Labial/Lingual/Buccal |

| | Order # | Head Size | Grit | Usage |
|---|------------|-----------|----------|--|
|  | W17DMHP-RA | 050 | M | Polisher - Labial/Lingual/Buccal |
|  | W17DHP-RA | 030 | F | Polisher - Labial/Lingual/Buccal |
|  | 8801LDF-FG | 014 | F | Diamond - Endo Access/ Small Adjustments |
|  | 8369DF-FG | 025 | F | Diamond - Occlusal Surface |
|  | 8856DF-FG | 016 | F | Diamond - Contact Surfaces/ Crown Removal |

TIPS & TECHNIQUES: Polishers produce best results when using all three steps in order (coarse, medium and fine); Use Copious amounts of water when using the diamonds.

Complimentary items: LD Grinders for heatless grinding (LD11C/M)

Power Suggestion: Use a Brasseler Forza FI Electric Attachment at low speed for polishers and Brasseler Forza F5 high speed attachment for diamonds.

Sterilization: Autoclave only

Recommended Speeds: 5,000-7,000 rpms

By Your Side In Dentistry



BRASSLER USA
DENTAL INSTRUMENTATION