



**BRASSELER**  
USA®

BY YOUR SIDE

 **enLITE**  
LED CURING LIGHT

- ▶ 360° Swivel Light Guide
- ▶ 420~490 nm Wavelength
- ▶ Digital Display
- ▶ Lightweight, Cordless
- ▶ Powerful, Lithium-Ion Battery



▶ Order online at [BrasselerUSA.com](http://BrasselerUSA.com) or call 800.841.4522.  
In Canada call 800.363.3838.



©2020 Brasseler USA. All rights reserved. Printed in the U.S.A.


B-5316-02.20

## ▶ 3 CURING MODES

### 1 | FULL POWER MODE

Full power for optimal high speed polymerization.


**CYCLE DURATION** Default 5 secs

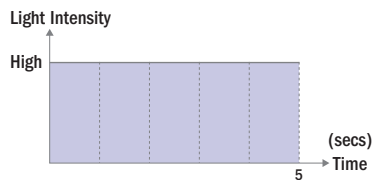

### 2 | RAMP-UP MODE

Light intensity gradually increases to full power to reduce composite shrinkage during polymerization.


**CYCLE DURATION** Default 10 secs

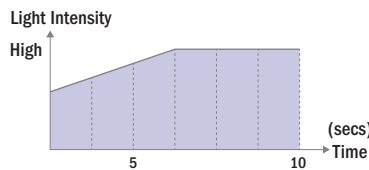

### 3 | PULSE MODE

Unit turns on and off every half second to reduce heat production.

**CYCLE DURATION** Default 15 secs

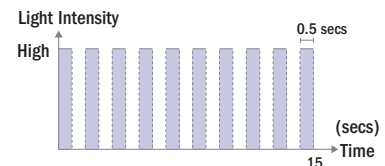


## ▶ ACCESSORIES

Order Number	Item Name
5027698U0	Light Guide*
5027969U0	Box of Barrier Sleeves
5028039U0	Anti-Glare Sleeve*
5028040U0	Anti-Glare Cone*
5027701U0	Power Adapter*

\*Included in complete set

## ▶ SPECIFICATIONS

<b>LED Lamp:</b>	5W
<b>Light Intensity:</b>	1,600 mW/cm
<b>Wavelength Range:</b>	420~490 nm
<b>Probe Size:</b>	8 mm
<b>Handpiece Weight:</b>	160 g
<b>Warranty:</b>	1-Year

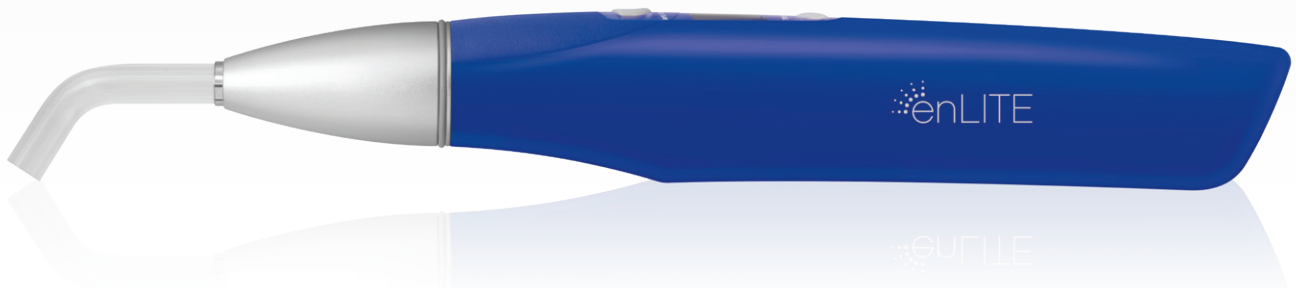

▶ Order online at [BrasselerUSA.com](http://BrasselerUSA.com) or call 800.841.4522.  
In Canada call 800.363.3838.

