



BRASSELER
USA®

BY YOUR SIDE

ComposiPro™ Brushes

One-Step Polishing System

Impregnated with silicon-carbide particles, these brushes create a highly polished surface without the need for messy polishing pastes. Their unique and versatile design makes it easy to polish difficult-to-access tooth surfaces and reduces chair time during your composite, porcelain, and metal restorations.

- ▶ One-step polishing system
 - ▶ Made of fine nylon fibers impregnated with silicone-carbide particles
 - ▶ Versatile; can be used on the surfaces of composite, porcelain, and metal restorations
 - ▶ Create a high gloss polish on textured surfaces
 - ▶ Remove discoloration during prophylaxis
 - ▶ Require no messy pastes
 - ▶ Autoclavable
- ▶ Order online at [BrasselerUSA.com](https://www.BrasselerUSA.com) or call **800.841.4522**.
In Canada call **800.363.3838**.



ComposiPro™ Brushes

Regular Cup

- Polishes large occlusal surfaces



Order No.	D209001
Head Size (1/10 mm)	060
Cutting Length (mm)	5.5
Grit	Medium/Fine
21 - RA - Right-Angle Latch (22 mm)	

Qty/Pack 5 & 20

Recommended speed: 5,000-7,000 rpm



Small Cup


- Polishes occlusal, facial, anterior, and cervical areas
- Used in confined or difficult to reach areas such as pedodontic cases, limited opening ability, and posterior teeth



Order No.	D209002
Head Size (1/10 mm)	040
Cutting Length (mm)	6.5
Grit	Medium/Fine
21 - RA - Right-Angle Latch (22 mm)	

Qty/Pack 5 & 20

Recommended speed: 5,000-7,000 rpm

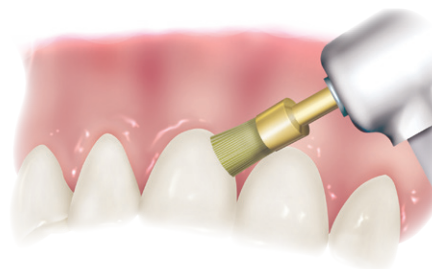

Point

- Polishes occlusal fissures and proximal areas



Order No.	D209003
Head Size (1/10 mm)	043
Cutting Length (mm)	6.5
Grit	Medium/Fine
21 - RA - Right-Angle Latch (22 mm)	

Qty/Pack 5 & 20

Recommended speed: 5,000-7,000 rpm



► Order online at BrasselerUSA.com or call **800.841.4522**.
In Canada call **800.363.3838**.

